

Under-Counter Pharmacy Refrigerator



HYC-118A

Scope of Application:

Versatile, quiet, under-counter pharmacy refrigerator features safety lock, digital display with multiple alarms. Widely used for storage of biological products, vaccines, drugs and reagents in within pharmacy, pharma manufacture, clinical, public health and other life science laboratories.

Innovative Design

- ◉ LED Display
- ◉ Ergonomic Design
- ◉ Multiple Alarms
- ◉ Lockable Design



Product Advantages



Superior Cabinet Temperature Uniformity, Auto Defrost

There are 9 test points within the interior which ensure superior temperature uniformity



Multiple Alarm

Capable of alerting failures due to high and low temperature, sensor error, power failure, low battery, door ajar and remote alarm



Energy Efficient Hydrocarbon (HC) Refrigeration System

Optimised refrigeration system, utilising high-efficiency HC compressor significantly reduces energy consumption



LED Display

Microprocessor control, digital display and temperature adjustment with an increment of 0.1°C



Ergonomic Design

810mm cabinet height makes it easy to place under test bench, saving storage space



Greener by Design


HC refrigerants are natural non-toxic and more environmentally friendly



Lockable Design


Safety lock to prevent unauthorized access

Product Introduction

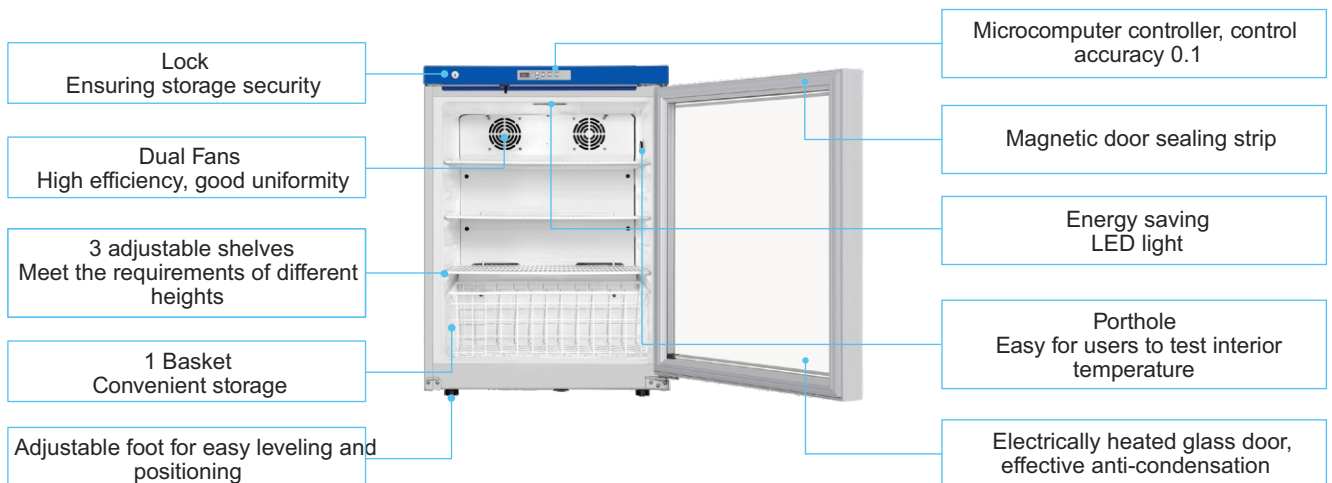


Temperature Curve

Cooling Down Curve

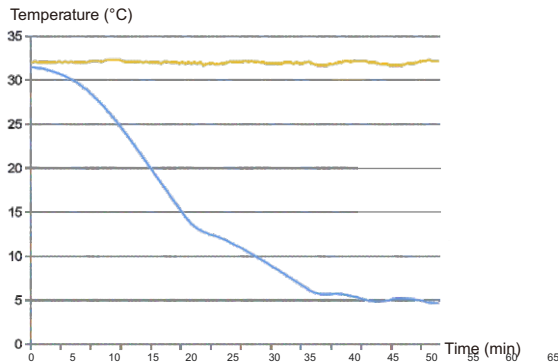

Stability Curve

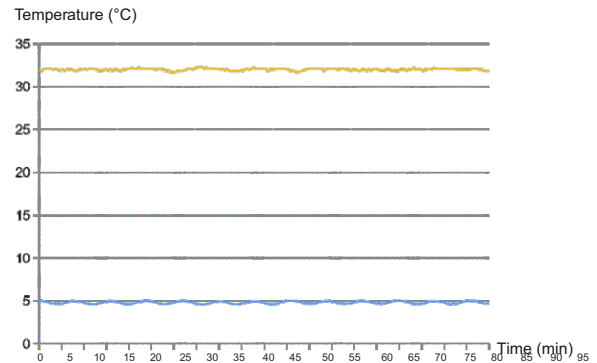

Specifications

— Internal Temperature
— Ambient Temperature

Model		HYC-118A		
Technical Data	Cabinet Type	Upright, glass door		
	Climate Class	N		
	Cooling Type	Forced air cooling		
	Defrost Mode	Auto		
	Refrigerant	HC		
	Noise (dB (A))	41		
Performance	Temp Range (°C)	2-8		
Control	Controller	Microprocessor		
	Display	LED		
Electrical Data	Power Supply (V/Hz)	220~240/50/60	115/60	
	Power (W)	230	230	
	Electrical Current (A)	1.5	1.5	
Construction	Capacity (L/Cu.Ft)	118/4.2		
	Net/Gross Weight (approx)	(kg)	46/51	
		(lbs)	101.5/112.5	
	Interior Dimensions (W*D*H)	(mm)	515*415*630	
		(in)	20.3*16.3*24.8	
	Exterior Dimensions (W*D*H)	(mm)	597*635*810	
(in)		23.5*25.0*31.9		
Packing Dimensions (W*D*H)	(mm)	680*690*910		
	(in)	26.8*27.2*35.8		
Accessories	Caster	Y		
	Foot	Y		
	Porthole	Y		
	Remote Alarm	Y		
	Shelves/Drawers	3/1		
Alarms	High/Low Temp	Y		
	Power Failure	Y		
	Sensor Error	Y		
	Low Battery	Y		
	Door Ajar	Y		
Certifications	CE	Y	/	
	UL	/	Y	

*Haier Biomedical reserves the right to change products and specifications without prior notice.

Marketed by:

BIOLINE TECHNOLOGIES

Sales Office : 404/405 Shanti dham heights, 90 feet Road
Parsik Nagar, Kalwa (W), Thane - 400 605, Maharashtra India,
Telefax : +91-22-2545 2111, +91-9594106999
Email : Info@bioline.in | **Web** : www.bioline.co.in
Branch : Delhi | Surat | Bangalore | Chennai | Trivandrum